



BleaseSirius

PRECISION. RELIABILITY. SIMPLICITY.



Developed and refined in partnership with you, Spacelabs Healthcare Ventilators offer the latest in patient-centered ventilation, addressing your real-life clinical needs. Our unique approach is based on proven critical care technology, enhancing ventilation across your range of patients and procedures.

You can choose the ventilator configuration and the way you interact with the Blease900 or Blease700 ventilators, both available with the BleaseSirius. A range of ventilation modes and waveform displays are available, enabling you to match your ventilator to your patient's changing requirements and your clinical and budgetary needs.

All Blease700 series ventilators offer Volume Control, Precision Pressure Control™ Ventilation (PPCV), and SIMV with Pressure Support. PPCV can enhance the performance, safety and precision of the ventilator by actively controlling flow into and out of the ventilator drive system throughout the entire respiratory cycle.

The Blease900 series offers all the features of the Blease700 series, plus Advanced Pressure Support™ that matches each patient's breathing



cycle on an 8.4" full color display. Ventilation monitoring is equally important to track changes or complications during ventilation. Measured at the patient utilizing Proximal Flow Sensing, the ventilator can display, store and compare pressure and flow volume loops.

Reduced Cost of Ownership

With over 80,000 units in use around the world, BleaseDatum™ vaporizers were the first to be free of routine maintenance requirements. We are so confident in these products that they are backed by a 10-year warranty* at no additional cost to you.

* 5-year warranty for Halothane.



Simple to complex patients, Spacelabs Healthcare offers you the choices you need, in terms of configuration and technique. We ensure them to satisfy your clinical preferences.

Delivering and ventilation is a long history of developing new technology to offer you the best in ventilation and precision, based on precision, accuracy and safety.

For you, including the latest technology to enhance ventilation precision and accuracy.

to move.


to your patients

Precision in an anesthesia ventilator to advanced mode SIMV-PC.

Ensuring flow and pressure as substantial advantages these parameters inside the Spacelabs Healthcare pioneered for the monitoring and control.

BleaseDatum vaporizers meet your requirements and are backed by a 5-year warranty.

Spacelabs Healthcare's innovative monitoring solutions provide vital signs monitors.

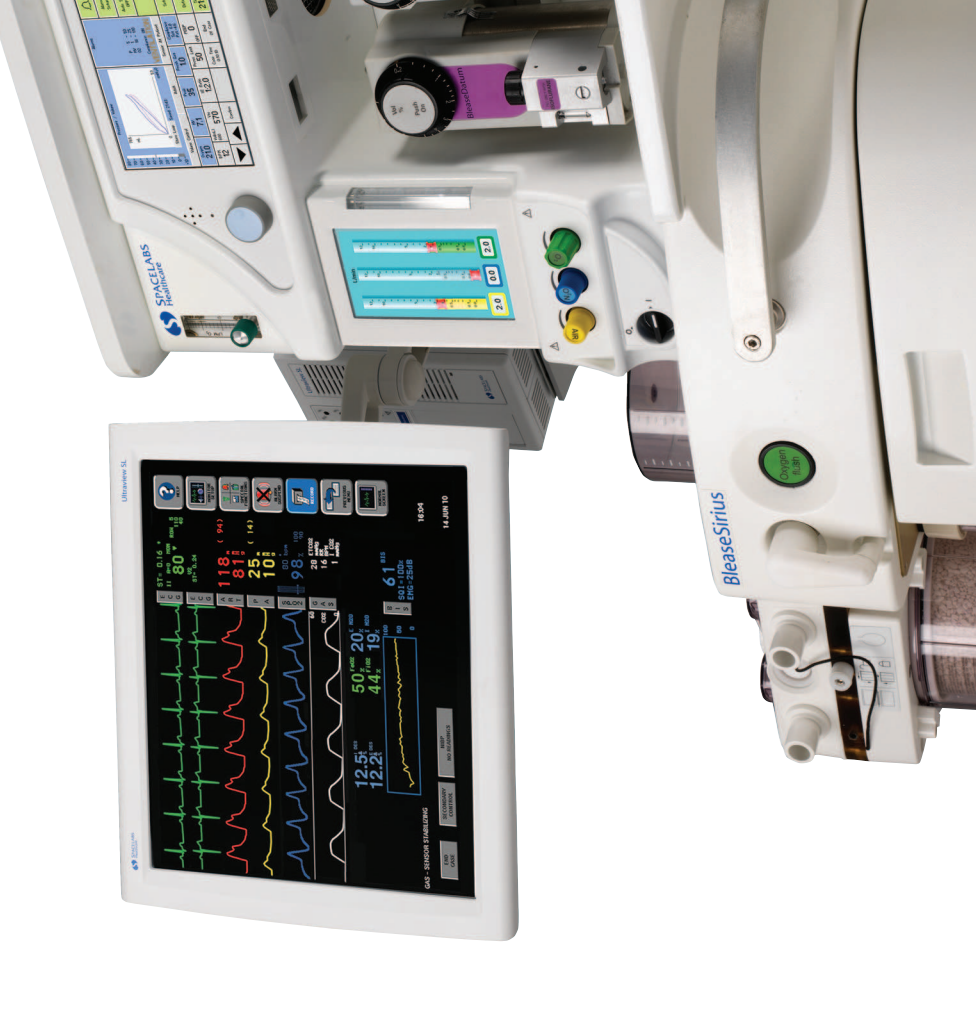

Electronic Flow Meter Option

- Instantly know the fresh gas flow to your patient
- Traditional flow control and positioning of information
- Individual gas data available for export

Identifying key information quickly and easily is critical to your practice.

The dedicated display is positioned where you would find traditional

flow tubes, allowing you to instantly know the fresh gas flow to your patient.

With traditional flow control for the Electronic Gas Flow Measurements



When selecting an anesthesia delivery company, consider **Spacelabs' 60 years of experience and our heritage of innovative and intuitive products** that have been developed in partnership with our customers to address real-life clinical needs. Our Customer First philosophy reflects our commitment to superior care for you and your patients and all of our products are supported by an award winning service team.

Adapt the equipment to your patients and procedures

Frame, ventilator and mounting options provide the customization capabilities you have asked for.

High flow or low flow, circle or semi-open systems, BleaseSirius adapts. An easily accessible scavenging connection point and the option of a highly visible, clearly-marked secondary common gas outlet, let you convert to a non-rebreathing circuit quickly and without contaminating the working environment.

The adjustable bag arm and APL valve, allow adjustment while maintaining focus on the patient.

BleaseSirius will integrate with your existing monitoring or, for the complete workstation experience, a range of innovative Spacelabs Healthcare anesthesia monitors.



www.spacelabshealthcare.com



Spacelabs Healthcare Limited
1 Harforde Court, John Tate Road
Hertford, SG13 7NW, United Kingdom

Spacelabs Healthcare Company Headquarters
5150 220th Ave. SE, Issaquah, WA 98029
PO Box 7018, Issaquah, WA 98027-7018

© 2010 Spacelabs Healthcare
Specifications subject to change without notice
A1001-oc10 / Drg No. 078-0477-00 Rev B

Telephone: +44 (0)1992 507700
Fax: +44 (0)1992 501213

Telephone: +1 425-657-7200
In North America: +1 800-522-7025
Fax: +1 425-657-7212